

Milikan Oil Drop Experiment

Jason Morgan

Department of Physics, Old Dominion University, Norfolk VA 23529

(Dated: April 29, 2014)

This experiment determines the charge of an electron by comparing the gravitational force on the oil drop to the electromagnetic force on the oil drop.

1. INTRODUCTION

In order to find the charge of the electron, the charge of several electron clusters were measured, then a common factor was calculated. The electron clusters were measured using a Millikan oil drop apparatus. This device suspends a small mist of charged oil droplets within a capacitor and allows an observer to time their movements with or without an electric field. By comparing the gravitational force to the electromagnetic force the complete charge on a single droplet can be measured. After several droplets are observed (and the smallest common denominator determined) the number can be compared to a known value.

2. EXPERIMENTAL SETUP AND PROCEDURES

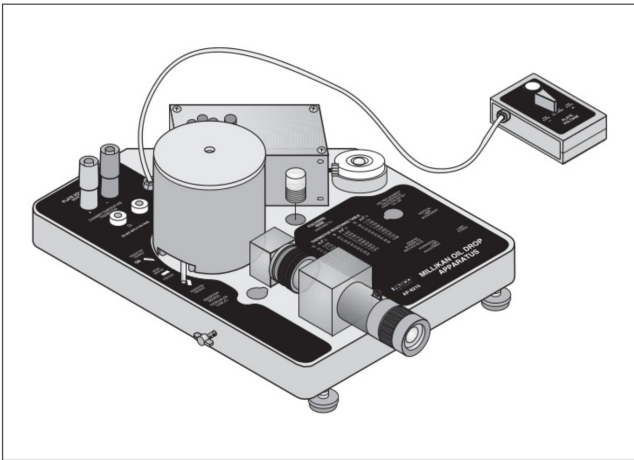


FIG. 1: The 'PASCO AP-8210' Millikan Oil Drop Apparatus

The device used in this experiment (the PASCO AP-8210) varies little from Millikan and Fletcher's original design. The chamber into which the oil droplets are sprayed consists of a cylindrical collection area with a small hole in the top just large enough for the tip of the atomizer. At the bottom of the chamber is a diverter which channels the oil mist to a minuscule viewing region. Attached to this region is a horizontal microscope lens, transparent grid, and light source. Above and below the viewing area are the parallel plates of the capacitor (which produce the electric field). Adjacent to the cham-

ber is an ionization source (which supplies the drops with electrons).

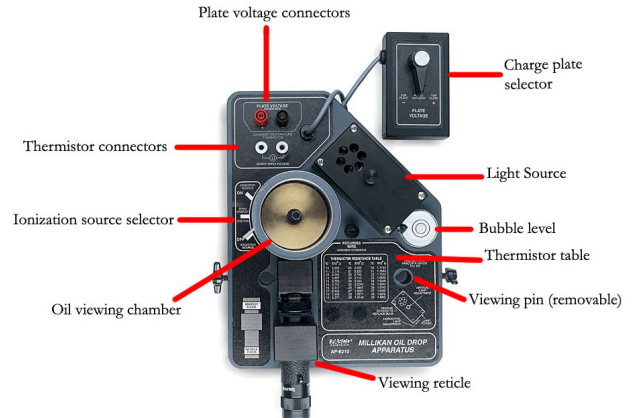


FIG. 2: Oil Drop Apparatus Components

This arrangement is mounted onto a flat base along with wiring for the plates and ionization source, a bull's eye level, and adjustable legs. The apparatus is also equipped with a thermistor in order to accurately gage the air temperature within the chamber (an important variable in air viscosity). Unlike the control for the ionization source, the control for the capacitor is attached via cord so that the frequent switching does not interfere with the drops. Likewise, the whole device is settled on a cushion of sand to minimize vibration transferring through the table.

To begin, the chamber was cleaned of oil, a small pin was paced in the viewing area and the light source was turned on. By adjusting the microscope and angle of light, the pin was brought into focus and removed. After returning the chamber to order, a continuous mist of oil drops were carefully introduced and the ionization source activated. Once droplets were visible from the viewing lens, the mist was discontinued and the ionization source turned off.

For a brief while the mist was observed to locate an optimal drop. Then the capacitor was alternately turned on and off. As this was done, the drop alternately rose and fell. Each transit was timed and recorded. Eventually, the drop would be lost to the background, measurements would cease and the process repeated. Each individual run of the experiment produced an independent set of data (along with the current values for the capacitor voltage and thermistor). The data was then

averaged and the charge for that particular drop was calculated using:

$$q = \frac{4}{3}\pi\rho g \left[\sqrt{\left(\frac{b}{2p}\right)^2 + \frac{9\eta v_f}{2g\rho}} - \frac{b}{2p} \right]^3 \frac{v_f + v_r}{E v_f}, \quad (1)$$

where:

q - total charge (*coulombs*)
 d - separation on plates (*m*)
 ρ - density of oil (*kg/m³*)
 g - acceleration of gravity (*m/s²*)
 η - viscosity of air (*Ns/m²*)
 b - barometric pressure (*Pa - m*)
 v_f - velocity of falling drop (*m/s*)
 v_r - velocity of rising drop (*m/s*)
 E - voltage/d ,

which concluded the experiment.

3. RESULTS AND DISCUSSION

Numerous attempts were made to accurately time the motion of a droplet, three were finally measured with some level of confidence (see Table I).

TABLE I:

Total Charge (coulombs)	Electrons
$1.31 \times 10^{-19} (\pm 4 \times 10^{-21})$	0.82
$3.02 \times 10^{-19} (\pm 2 \times 10^{-21})$	1.88
$6.72 \times 10^{-19} (\pm 1 \times 10^{-21})$	4.19

The smallest charge of the three (which appeared to possess only one excess electron) held a charge of $1.31 \times 10^{-19} (\pm 4 \times 10^{-21})$ coulombs, which varies from the expected value of 1.60×10^{-19} coulombs by 22%.

The percent difference from expected values predictably decreased with both the second and third data sets, as the number of electrons increased, giving around 5% differences. Without knowledge of the accepted value for an electron charge, the third data set would have been calculated as possessing around 5 electrons instead of 4, and the second set would have been calculated at closer to 2.5 electrons (which starkly contradicts quantization).

4. CONCLUSION

This was a great experiment to perform. While it was quite difficult to achieve good results it was a good demonstration of one can measure the charge of a charged particle.

-
- [1] Sunny Bishop, *MILLIKAN OIL DROP APPARATUS*, (10101 Foothills Blvd. Roseville, CA).
 [2] Wikipedia. "Oil drop experiment.", *Wikipedia*, Wikimedia Foundation, 21 April 2014. Web. 22 April 2014.

- [3] Wikipedia. "Robert Andrews Millikan.", *Wikipedia*, Wikimedia Foundation, 20 April 2014 . Web. 22 April 2014.